P0#1SB004143_GBS-2232-F_PN15313_说明书

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QUICK START GUIDE

Terraillon®

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

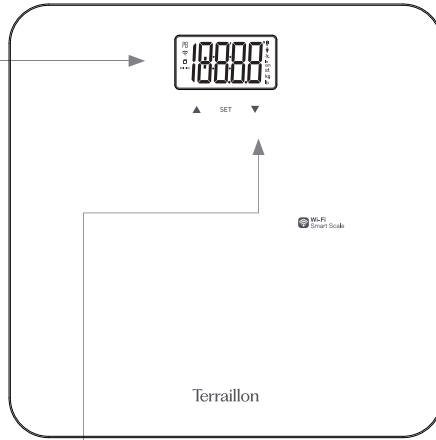
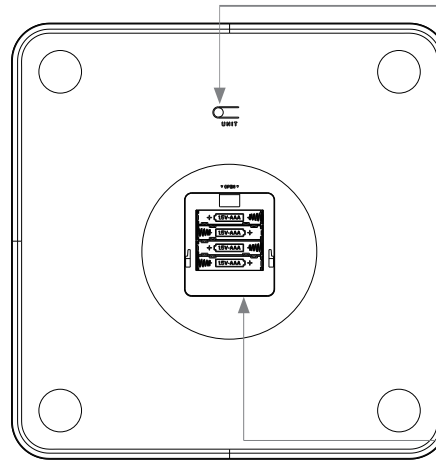
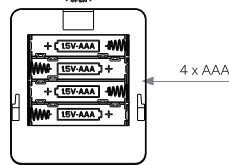
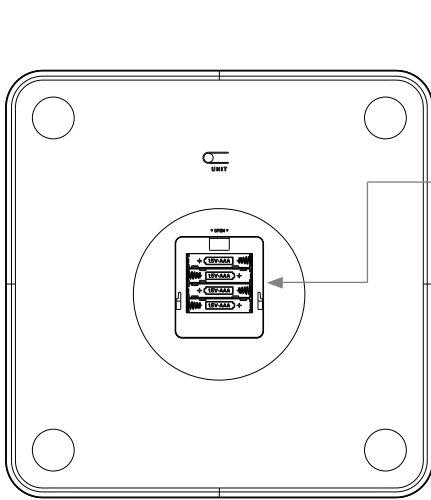
INSERT BATTERIES INTO THE DEDICATED COMPARTMENT.

- | | |
|--|---|
| FR Insérez les piles dans le compartiment dédié. | EL ΤΟΠΟΘΕΤΗΣΤΕ ΤΙΣ ΜΠΑΤΑΡΙΕΣ ΣΤΟ ΕΙΔΙΚΟ ΤΜΗΜΑ. |
| NL Plaats de batterijen in het daarvoor bestemde vak. | SV SÄTT I BATTERIERNA I DET AVSEDDA FACKET. |
| ES Inserte las pilas en el compartimento específico. | NO SETT BATTERIERNE INN I DET DEDIKERTE ROMMET. |
| IT Inserire le batterie nel vano apposito. | FI ASETA PARISTOT NIILLE TARKOITETTUUN LOKEROON. |
| PT Insira as pilhas no compartimento dedicado. | DA SÆT BATTERIERNE I DET DERTIL BEREGNEDE RUM. |
| DE Legen sie die batterien in das entsprechende fach ein. | |

GO OVER TO YOUR INTERNET ROUTER AND MAKE SURE YOU KNOW YOUR Wi-Fi PASSWORD. THEN CONNECT YOUR SMARTPHONE TO YOUR Wi-Fi NETWORK.

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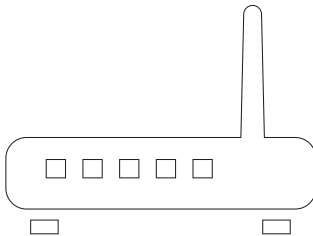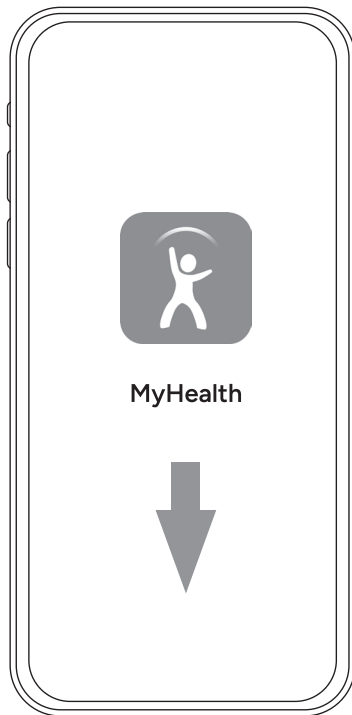DOWNLOAD MyHealth MOBILE APPLICATION.

- FR** Téléchargez l'application mobile MyHealth.
- NL** Download de mobiele applicatie MyHealth.
- ES** Descarga la aplicación móvil MyHealth.
- IT** Scaricare l'app mobile MyHealth.
- PT** Descarregue a aplicação móvel MyHealth.
- DE** Laden Sie die mobile App MyHealth herunter.
- EL** ΚΑΤΕΒΑΣΤΕ ΤΟΝ ΠΡΟΠΟΝΗΤΗ ΕΥΕΞΙΑΣ – ΤΗΝ ΕΦΑΡΜΟΓΗ MY HEALTH MOBILE.
- SV** LADDA NED WELLNESS COACH – MY HEALTH MOBIL-APP.
- NO** LAST NED MOBILAPPEN WELLNESS COACH – MY HEALTH.
- FI** LATAA WELLNESS COACH – MY HEALTH MOBILE -SOVELLUS.
- DA** DOWNLOAD WELLNESS-COACH – MY HEALTH-MOBILAPP.

CREATE YOUR ACCOUNT OR LOG IN WITH YOUR CREDENTIALS.

- FR** Créez votre compte ou connectez-vous avec vos identifiants.
- NL** Maak een account aan of meld u aan met uw inloggegevens.
- ES** Cree una cuenta o conéctese con sus credenciales.
- IT** Creare un account o accedere con i propri dati di accesso.
- PT** Crie a sua cunta ou ligue-se com os seus identificadores.
- DE** Erstellen Sie Ihr Konto oder melden Sie sich mit Ihrem Login an.
- EL** ΔΗΜΙΟΥΡΓΗΣΤΕ ΤΟ ΛΟΓΑΡΙΑΣΜΟ ΣΑΣ Ή ΣΥΝΔΕΘΕΙΤΕ ΜΕ ΤΑ ΔΙΑΠΙΣΤΕΥΤΗΡΙΑ ΣΑΣ.
- SV** SKAPA ETT KONTO ELLER LOGGA IN MED DIN INLOGGNINGSUPPGIFTER.
- NO** OPPRETT EN KONTO ELLER LOGG INN MED BRUKERNAVN OG PASSORD.
- FI** LUO TILI TAI KIRJAUDU SISÄÄN TUNNUKSILLASI.
- DA** OPRET DIN KONTO, ELLER LOG IND MED DINE BRUGEROPLYSNINGER.


QUICK START

FR Démarrage rapide - **NL** Snelstarten - **ES** Inicio rapido - **EL** ΓΡΗΓΟΡΗ ΕΚΚΙΝΗΣΗ - **SV** SNABBSTART - **NO** HURTIGSTART
IT Avvio rapido - **PT** Início rápido - **DE** Schnellstart - **FI** PIKA-ALOITUS - **DA** QUICKSTART

DURING THE CREATION OF YOUR ACCOUNT, AND WHEN THE APPLICATION INVITES YOU TO DO SO, SELECT THE WI-FI BATHROOM SCALES TO START THE PAIRING PROCEDURE.

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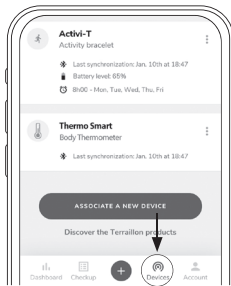
NOTES: IF YOU ARE CONNECTED WITH A PRE-EXISTING ACCOUNT, PRESS THE BUTTON AS INDICATED, THEN "ASSOCIATE A NEW DEVICE". WHEN THE APPLICATION INVITES YOU TO DO SO, SELECT THE WI-FI BATHROOM SCALES TO START THE PAIRING PROCEDURE.

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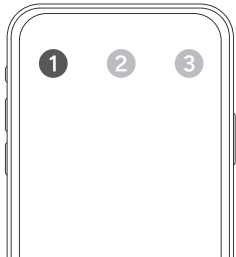FOLLOW THE INSTRUCTIONS ON THE MOBILE APPLICATION.

- FR** Suivez les instructions dans l'application mobile.
- NL** Volg de instructies in de mobiele applicatie.
- ES** Siga las instrucciones de la aplicación móvil.
- IT** Seguire le istruzioni dell'applicazione mobile.
- PT** Siga as instruções da aplicação móvel.
- DE** Folgen sie den anweisungen in der mobilen anwendung.
- EL** ΑΚΟΛΟΥΘΗΣΤΕ ΤΙΣ ΟΔΗΓΙΕΣ ΣΤΗΝ ΕΦΑΡΜΟΓΗ ΓΙΑ ΚΙΝΗΤΑ ΤΗΛΕΦΩΝΑ.
- SV** FÖLJ INSTRUKTIONERNA I MOBILAPPEN.
- NO** FØLG INSTRUKSJONENE PÅ MOBILAPPEN.
- FI** NOUDATA MOBIILISOVELLUKSEN OHJEITA.
- DA** FØLG VEJLEDNINGEN I MOBILAPPEN.

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EN

Please read this section carefully to familiarize yourself with features and operations before using the unit. The warning signs and the sample icons shown here are listed in order for you to use this product safely and correctly as well as to prevent product damage, risk and injury to you or others. This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.

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Wi-Fi 802.11 b/g/n (2.4 Ghz compatible) - 16.51dBm - 17.94dBm

Weight & BMI (5 kg-180 kg)

8 user memories (with automatic recognition)

Battery powered (4 x AAA alkaline batteries (not included))

Dimensions: 32,7 x 32,7 x 2,19cm

MAINTENANCE

Do not disassemble the scale: it contains no user serviceable parts. Do use a wet soft cloth, dipped into water and wrung out, to wipe the dirt. Then use a dry soft cloth to dry up the device. DON'T use propellant, abrasive or other chemicals.

PROTECT THE ENVIRONMENT

Dispose of used batteries in a designed container so that they can be collected and recycled. After life, dispose of this product at a designated waste recycling point.

WARNING

Do not expose batteries to high temperatures on account of the risk of liquefaction, vaporisation and explosion. Please ensure the batteries are replaced with batteries of the same type (otherwise your device may be damaged).

DECLARATION OF CONFORMITY

Hereby, Terrailon SAS declares that MASTER FORM ULTRA complies with the essential health and safety conditions. This product meets the essential requirements and other relevant provisions of the directive RED 2014/53/UE. The complete declaration of conformity is available on <https://terrailonhelp.zendesk.com/hc/en-us> (Declaration of Conformity section).

WARRANTY

This product is guaranteed against defects in materials or manufacturing for 2 years. During this period, any such defects will be repaired free of charge (proof of purchase must be shown in the case of a complaint under warranty). The warranty does not cover damage resulting from improper use or negligence. In the event of a claim, first contact the shop where you purchased your product.

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EN If you have any difficulties or need more information, you can find our help centre by scanning the QR code opposite with your smartphone. If you cannot find the answer to your question, you can also contact our After Sales Service: support.uk@terraillon.com.

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EN

FR

Terraillon SAS France & Headquarters

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Terraillon Asia Pacific Ltd

4/F, Eastern Centre
1065 King's Road
Quarry Bay – Hong Kong
Tel: +852 (0)2960 7200
customerservice@terraillon-asia.com

Terraillon Corp USA

contact@terraillon.fr

Terraillon UK

support.uk@terraillon.com

Terraillon® CE

MASTER FORM **ULTRA**