WEBCOACH PRIME

USER MANUAL | GUIDE D'UTILISATION

NL Handleiding | IT Manuale d'uso | ES Manual del usuario

DE Bedienungsanleitung | PT Manual do utilizado

Terrailon®

www.terrailon.com

Made for / Compatibilité


iPhone® 4S/5/5C/5S/6/6Plus
iPod® Touch 5th generation

iPad® 3/4/Air
iPad® Mini

Android 4.3
Bluetooth Smart 4.0 / Bluetooth Smart Ready

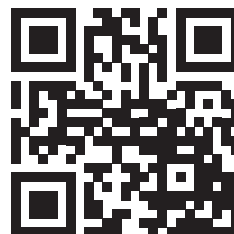iPhone is registered trademark of Apple Inc.

Android is a trademark of Google Inc.

The Bluetooth® word mark and logos are registered trademarks
owned by Bluetooth SIG, Inc.

Other trademarks and trade names are those of their respective owners.

WEBCOACH PRIME

EN	TERRAILLON CONNECTED BODY COMPOSITION ANALYZER INSTRUCTION MANUAL	4
FR	IMPÉDANCEMÈTRE CONNECTÉ TERRAILLON NOTICE D'UTILISATION	8
NL	TERRAILLON LICHAAMSVETWEEGSCHAAL HANDLEIDING	12
IT	BILANCIA MASSA GRASSA TERRAILLON MANUALE DI ISTRUZIONI	16
ES	BÁSCULA DE MASA GRASA DE TERRAILLON MANUAL DE INSTRUCCIONES	20
DE	TERRAILLON KÖRPERFETTANALYSEWAAGE BEDIENUNGSANLEITUNG	24
PT	BALANÇA DE MASSA GORDA TERRAILLON MANUAL DE INSTRUÇÕES	28

Thank you very much for selecting the Terrailon Glass Body Fat Analyzer
WEB COACH PRIME.



Don't use this product if you wear a pacemaker or implantable electronic device • Any information provided by this device is in no way meant to treat, cure or prevent any disease or illness from happening.

FIND MORE INFORMATION ON

<http://www.terraillon.com>

PRODUCT OVERVIEW

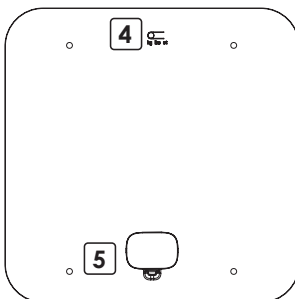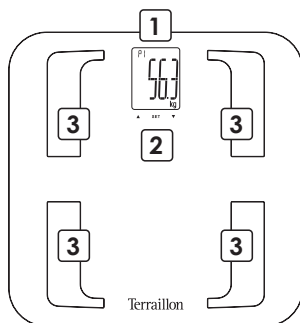- 1** LED Display
- 2** Profile selection
- 3** Electrodes
- 4** Unit and pairing button
- 5** Battery compartment


FIRST USE

- A** Insert the batteries (3xAAA)
- B** Download the Terrailon Coach application


Download the Terrailon
Wellness Coach
user manual to know more
about our app


C Create an account

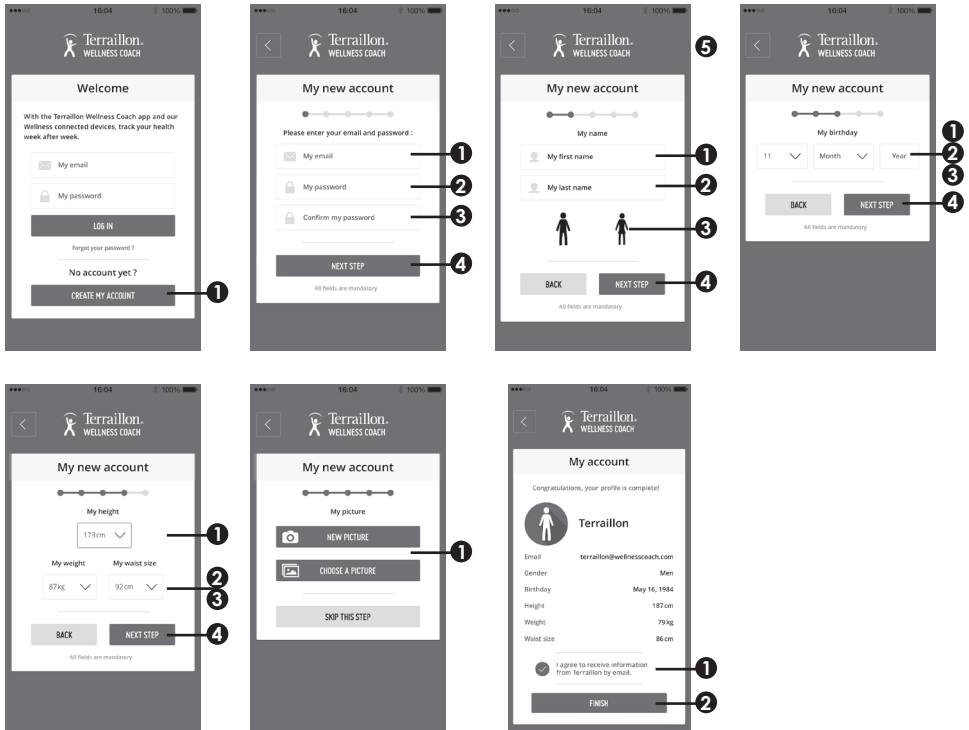D Turn on Bluetooth on your Smartphone (Setting > Bluetooth > ON)

E Hold "UNIT" button in the back to start pair-up

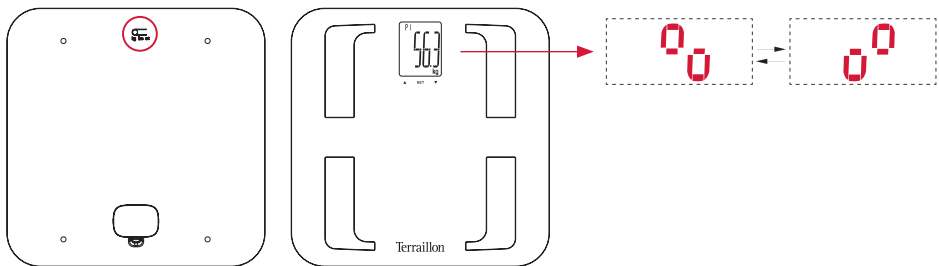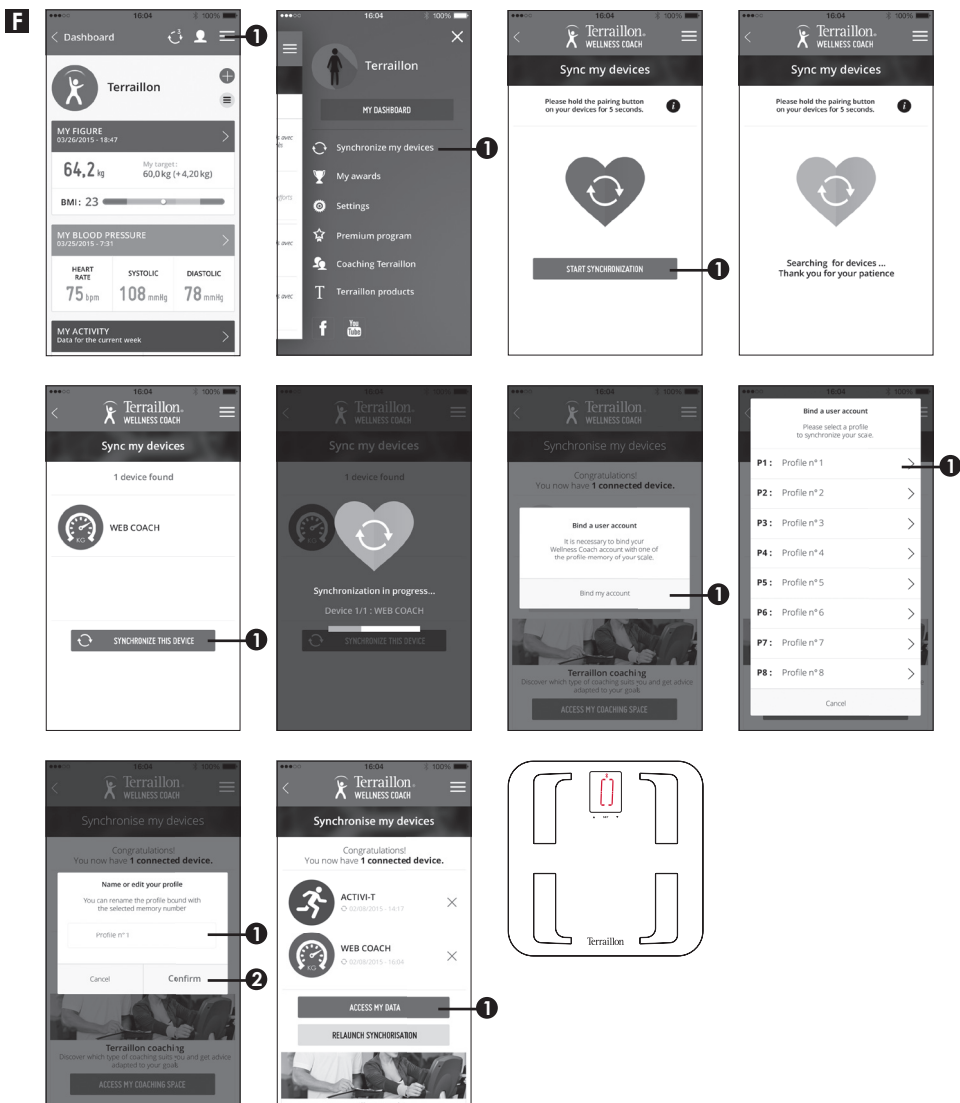- G** It's important to bind your account with the scale for automatic recognition of user during weighing. You can pair the scale with 8 smartphones maximum. To do that, just repeat steps B to F on an other smartphone and select another profile number. To bind an other user on the same smartphone, please log in with the other user account and repeat steps E and F. Then select another profile number. Before you first weighing, please select your profile number on the scale.

If two profiles have a similar weight, the scale will ask you to confirm your profile number. On the scale, press the left arrow ▲ if your profile number is the one on the left or press the right arrow ▼ if your profile number is the one on the right. On the below example, if you are profile 3, you have to press the left arrow ▲.

START MEASURING

- A** The scale will automatically switch on as you step on the platform barefooted. The symbol ✱ remains fixed on the screen, that means that the connection with the app is done. If it blinks, please be sure that, on your Smartphone, your Wellness Coach app is open and the Bluetooth is on. Stand still and keep full contact with the electrodes until the LED stops blinking "BIA". The symbol ☒ blinks on the screen and disappears. It means that the data has been transmitted to the app. If the symbol remains on the screen, it means that data transmission has failed. The data is saved in the scale. Please check Bluetooth connection (symbol ✱).

Be sure to always weigh yourself on a hard and flat surface.

MAINTENANCE

Use water only. DON'T use propellant or other chemicals.

WARRANTY & ENVIRONMENT PROTECTION

This product is guaranteed for 2 years against material and manufacturing defects. During this period, such defects will be repaired free of charge (proof of purchase must be presented if the scale is under warranty). This warranty does not cover damage resulting from accidents, improper use or negligence. If you have a complaint please contact the store where you purchased the scale first.

Waste electrical products should not be disposed of with household waste. Please recycle where facilities exist. Check with your local authority or retailer for recycling advice.

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Terrailon SAS France & Headquarters
Terrailon UK Ltd
Terrailon Asia Pacific Ltd
Terrailon Corp USA