Petit-Poids

FR	Pèse-bébé électronique Terraillon	p.02
EN	Terraillon electronic baby scale	p.04
DE	Terraillon elektronische Babywaage	p.06
IT	Pesabambini elettronica Terraillon	p.08
ES	Pesabebés electrónico de Terraillon	p.10
NL	Terraillon elektronische babyweegschaal	p.12
PT	Balança de bebé electrónica Terraillon	p.14

Terraillon®

CE

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
*Dear Customer, thank you for purchasing this **electronic Terraillon baby scale**.
We hope you enjoy using the
product. To ensure your complete satisfaction, we recommend
you read these instructions carefully before using it.*

Terraillon electronic baby scale

This Terraillon electronic baby scale has been designed based on high quality and safety standards. Its **measurement freeze** function means that the weight is displayed, even if the baby keeps moving. Its highly ergonomic design guarantees comfort and stability, and it can memorise a measurement and calculate the difference between the current and previous weights.

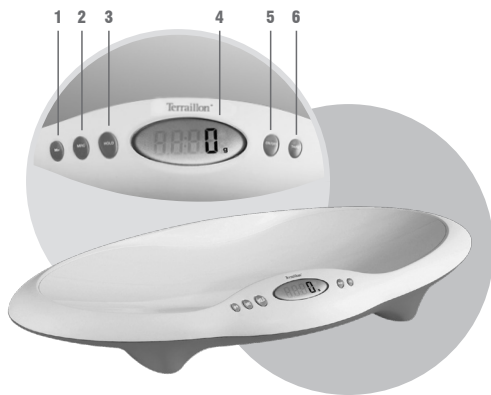SPECIFICATIONS

- 1- Memorisation of weight (M+)
- 2- Recall of memorised weight (MRC)
- 3- Weight freeze function (HOLD)
- 4- Large 78.5 x 37 mm LCD screen
- 5- On/Off switch (ON/OFF)
- 6- Automatic tare (TARE)

FEATURES

- Maximum weight: 15 kg
- Increments: 10 g
- Measurement freeze function
- Weight memory function
- Indicates the difference between the current and previous weights
- Automatic tare
- Power supply: Two 1.5V AAA batteries, LR03 type (included)
- Low battery indicator

OPERATION

1- PUTTING IN THE BATTERY

The appliance uses two 1.5V AAA batteries (LR03 type). Insert the two 1.5V batteries into the battery compartment underneath the baby scale. Respect the polarities +/-.

2- USE

Press the ON/OFF button. All the elements on the screen are displayed for a few seconds, followed by "0.00kg" and the symbol "⌚". When the symbol "⌚" disappears and only "0.00kg" appears on the screen, place the baby in the cradle.

If your baby keeps still, his/her weight will automatically stabilise after a few seconds and then appear on the screen.

If your baby is moving, the symbol "⌚" appears on the left-hand side of the screen and a non-stabilised weight scrolls on the screen. Press the HOLD button to freeze the baby's weight. The "⌚" symbol is replaced by "HOLD".

The screen displays the baby's stabilised weight.

When it is not being used, the baby scale automatically switches off after approximately 20 seconds.

a- Automatic Tare Function

The automatic tare function allows you to weigh your baby comfortably lying on a towel without including the weight of the towel.

- Switch on the scale by pressing the ON/OFF button.
- Put the towel on the cradle.
- Press the TARE button, the screen displays "---" then "0.00kg".
- Place your baby comfortably on the towel.

Only the baby's weight is displayed.

(When you remove both the baby and the towel, a negative weight may appear. This is the corresponding weight of the towel in negative.)

b- Weight memorisation

When the weight of the baby appears on the screen, you can store it in the memory by pressing the "M+" button, and it will continue to be displayed for a few seconds. The "M+" symbol appears.

c- Recalling the recorded weight

Press the **ON/OFF** button to switch on the baby scale.

By pressing the "MRC" button, the previous memorised weight will appear on the screen with the sign "⊕" to the right. The memorised weight will continue to be displayed for a few seconds.

The recorded weight will remain in the memory until a new weight is recorded and whilst the battery continues working.

d- Difference between two weight measurements

The baby scale automatically displays the difference between the current weight and the previously memorised weight.

The difference in weight is displayed using the ▲ or ▼ signs, depending on whether the weight is higher or lower.

PROBLEMS	REASON	SOLUTION
Err2	The baby scale is overloaded (weight over 15 kg).	Remove your child immediately or you may damage the baby scale.
ErrE	Baby scale initialisation problem.	Switch the baby scale on and off with no weight on the scale.
bAtt	The battery is worn out.	Replace the battery immediately.

DIRECTIONS FOR USE

- Clean your baby scale with a damp cloth without letting water get into the appliance.
- Remove the two 1.5V batteries if the baby scale is not going to be used for lengthy periods..

Used batteries should be placed in the appropriate collection bin so they can be recycled. Do not mix different types of batteries, or use new and old batteries at the same time. At the end of its serviceable life, this product must be disposed of at a recycling centre for electrical and electronic equipment.

WARRANTY

Your Terraillon electronic baby scale is guaranteed for two years against material and manufacturing defects.

During this period, such defects will be repaired free of charge (proof of purchase must be presented if the claim is made under warranty).

This warranty does not cover damage resulting from accidents, improper use or negligence.

In the event of a claim, please contact the store where you first bought the item.

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Terraillon SAS France & Headquarters
Terraillon UK Ltd
Terraillon Asia Pacific Ltd
Terraillon Corp USA