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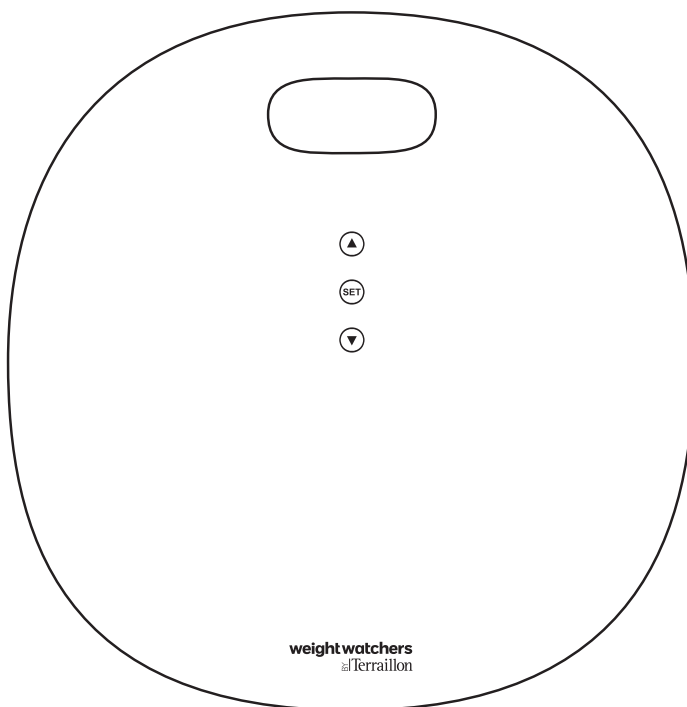
AUTOMATIC ELECTRONIC BATHROOM SCALE - INSTRUCTION MANUAL

Dear Customer,

Thank you for purchasing this product. We hope you will find it satisfactory. Please read these instructions carefully to take full advantage of your product.

WHAT IS EASY FOLLOW?

The EASY FOLLOW scale indicates:

- Your weight
- Your weight gain or loss compared to your previously saved weight
- Your Body Mass Index (BMI)

WHAT IS BMI (BODY MASS INDEX)?

• **For adults from 19 to 99 years*:**

It is the ratio between a person's weight and height. This index is used to assess a person's body fat. Calculation:

$$\text{BMI} = \text{weight (kg)} / [\text{height (m)}]^2$$

* For adults over 65, the BMI classification can sometimes be inaccurate. However, using the Body Control does not in any way represent a health risk.

> 30	Between 25 and 29,9	Between 18,5 and 24,9	< 18,5
Obese	Overweight	Normal	Underweight


READ BEFORE USE

- Do not disassemble the scale as incorrect handling may cause injury.
- The product is not intended for professional usage in hospitals or medical facilities; it is intended for home/consumer use only.
- Clean after use with a lightly dampened cloth without allowing any water to get into the device. Do not use solvents or immerse the product in water.
- Avoid excessive impact or vibration to the unit, such as dropping the product onto the floor.
- The scale must be reset each time it has been moved.
- Pause after each weighing.

AUTOMATIC ELECTRONIC BATHROOM SCALE - INSTRUCTION MANUAL

PRODUCT DESCRIPTION

WEIGHING PLATFORM

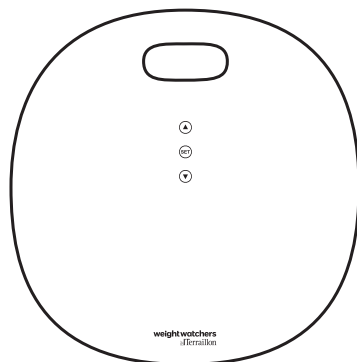FUNCTION KEYS

Up


Set button

Down

LCD DISPLAY

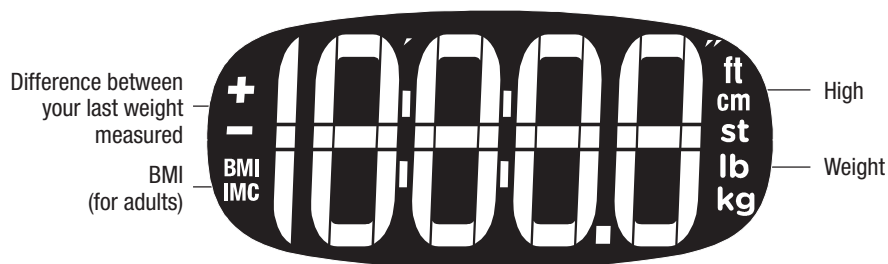weight watchers BY Terrillon®

AUTOMATIC ELECTRONIC BATHROOM SCALE - INSTRUCTION MANUAL

PREPARATION BEFORE USE

Insert the batteries

Open the battery cover on the back of the scale and the one on the back of the control terminal. Insert the batteries in both compartments. Be sure that the polarity of the battery is set properly for proper function (the polarity as indicated inside the battery compartment cover). If you do not intend to use this unit for a long period of time, it is advisable to remove the batteries before storage.

Initialisation

The scale needs to be initialised at the first time of use, each time it has been moved or after battery replacement. Press briefly and sharply with your foot in the middle of the platform. The display reads « 0.0 » then the scale automatically turns off. Your scale is now ready for use.

If your scale has a conversion switch, you can change the measuring units from metric to imperial or vice versa (i.e KG/ST/LB) by simply pressing the conversion switch located on the under-side of the scale.

INSTRUCTION MANUAL

When weighing yourself for the first time

Set up personal data (stature) before using, it can set up 4 users' data.

- Hold 3s "set" to enter the mode of setup:
 - The flashing Pn (n can be 1~4) means the memory of the data.
 - Press ▲ & ▼ to change number then press "SET" to confirm.
- Set up the stature: Press UP, DOWN, and then press SET to confirm.
- The scale will back to the mode of weighing after finishing setting up.
- Stand on the scales.
- The scales will indicate your weight as well as your BMI.

When weighing yourself subsequently

- Stand on the scales.

Auto Recognition

If user does not choose a profile (Pn) manually and step on the scale, the weight value will blink twice and be locked. After that, the product will auto-identify user. The product compares all users' last weighing measurement within +-3kg.

- If the product fails to identify then it will only show weight value.
- If product recognizes one user, it will directly indicate your current weight plus the difference between your last weight measured (a «-» sign in front of the figure indicates weight loss), and then your BMI.
- When it recognizes more than one user profile, it will show all it recognized. For example, the product recognizes 2 users, it will show the 2 user's number for 4 seconds, repeat 2 times. When the LCD blinks, user can press SET to confirm user no. and then the scale will indicate your current weight plus the difference between your last weight measured (a «-» sign in front of the figure indicates weight loss), and then your BMI.
- If user does not choose user profile during recognition, then the scale automatically turns off.

AUTOMATIC ELECTRONIC BATHROOM SCALE - INSTRUCTION MANUAL

Quick memory recall

- With the device switched off, press the [▲] or [▼] buttons to go to the desired memory slot and the SET to confirm your choice.

WARNING INDICATORS

Incorrect initialisation - The scale must be reset.

Overload: remove weight immediately.

Batteries are low: Replace batteries

PROTECT THE ENVIRONMENT

Dispose of used batteries in a designated container so that they can be collected and recycled. Do not mix different types of batteries. Do not mix new and used batteries.

After life, dispose of this product at a designated waste recycling point.

GUARANTEE

This product is guaranteed against defects in materials or manufacturing. During the guarantee period, any such defects will be repaired free of charge (proof of purchase must be shown in the case of a complaint under guarantee). The guarantee does not cover damage resulting from improper use or negligence. In the event of a claim, first contact the shop where you purchased the product.

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© 2016 Weight Watchers International, Inc. Tous droits réservés.
Terraillon SAS France & Headquarters
Terraillon UK Ltd
Terraillon Asia Pacific Ltd
Terraillon Corp USA