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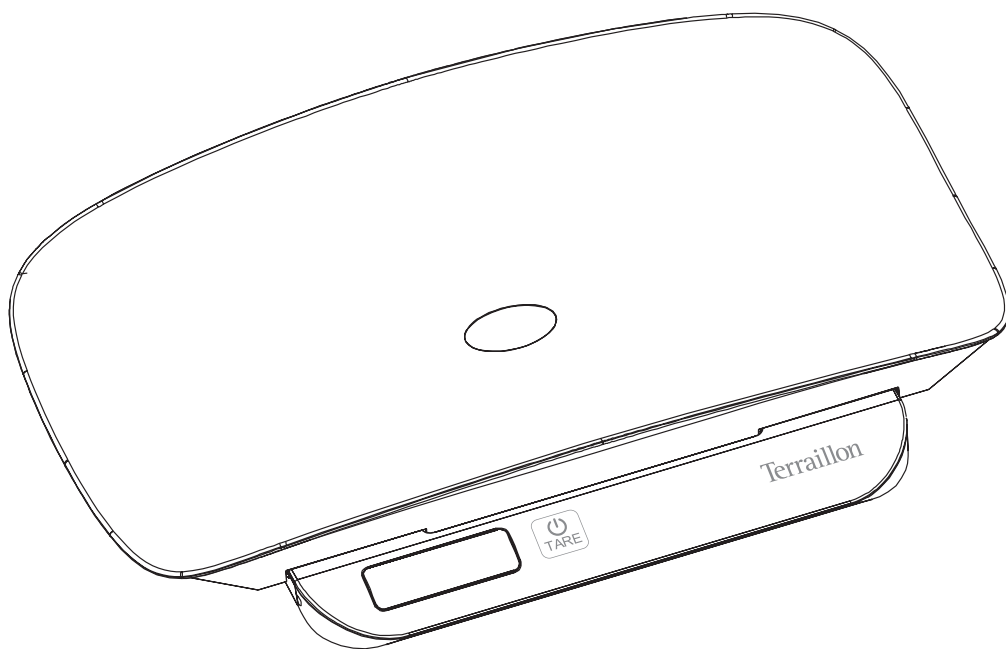
TERRAILLON EVOLUTIVE BABY SCALE


Dear Customer,

Thank you for purchasing this evolutive Terraillon baby scale. We hope you enjoy using the product. To ensure your complete satisfaction, we recommend you read these instructions carefully before using it.

READ BEFORE USE

- This scale is designed for weighting Baby and Toddler (up to 20 kg).
- Place the scale on a hard and flat surface.
- Do not disassemble the scale as incorrect handling may cause injury.
- The product is not intended for professional usage in hospitals or medical facilities; it is intended for home/consumer use only.
- Clean after use with a lightly dampened cloth without allowing any water to get into the device. Do not use solvents or immerse the product in water.
- Avoid excessive impact or vibration to the unit, such as dropping the product onto the floor.
- Do not press or put any object on platform when not use.
- Extreme electromagnetic influences can adversely affect the performance of the scale. Keep the scale at a sufficient distance from such equipment or turn off the equipment when using the scale. Once the scale is interfered, shut it off and turn it on again.
- Remove the batteries if the baby scale is not going to be used for lengthy periods.

DESCRIPTION OF THE PRODUCT

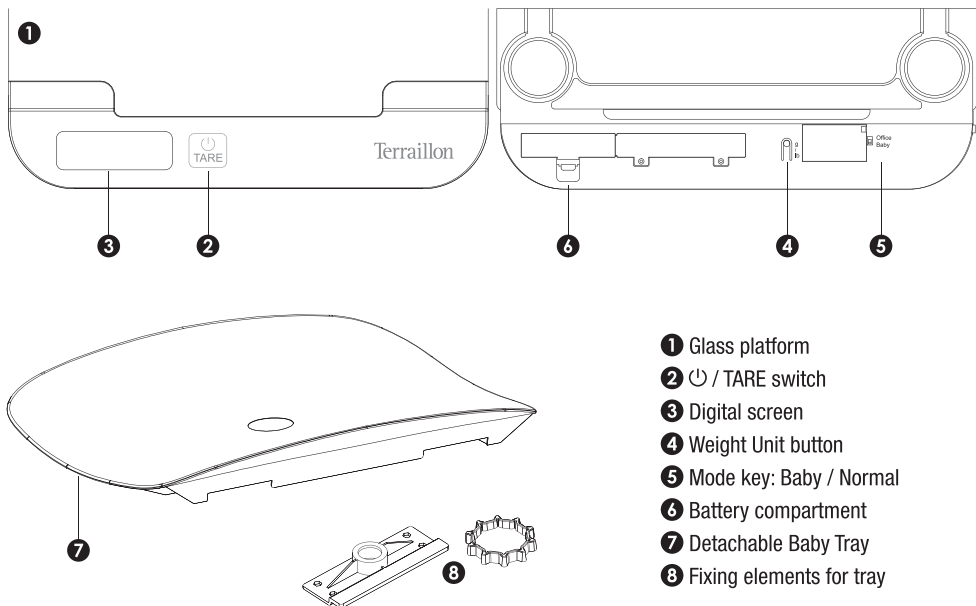- ① Glass platform
- ② ⏻ / TARE switch
- ③ Digital screen
- ④ Weight Unit button
- ⑤ Mode key: Baby / Normal
- ⑥ Battery compartment
- ⑦ Detachable Baby Tray
- ⑧ Fixing elements for tray

TERRAILLON EVOLUTIVE BABY SCALE

PRODUCT SPECIFICATION

Maximum weight: 20 kg / 44 lb /

Display: < 10 kgs, "G" display ; > = 10 kgs, "KG" display

Increments: 5 g / 0.2 oz

Normal Mode

Baby Mode with auto Weight Hold

Power supply: 3 x 1.5 V AA batteries (not included)

Low battery indicator

PREPARATION BEFORE USE


The appliance uses 3 X 1.5 V AA batteries (LR6 type). Insert the batteries into the battery compartment underneath the baby scale. Respect the polarities +/-.

Select the weight unit -g/lb/st (weight unit button on the rear side of LCD panel).

OPERATION

BABY MODE

1. Select first BABY MODE with the mode switch on rear side of LCD panel.
If you choose baby mode, "baby" text will show on LCD panel top right hand side.
2. Please put scale tray on the scale if you choose baby mode.

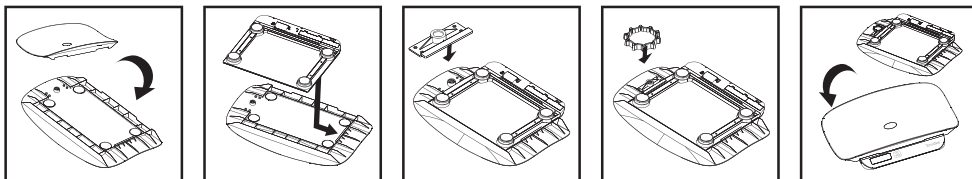Place the scale on a hard and flat surface.

3. Press the  / TARE button. All the elements on the screen are displayed for a 5 seconds, followed by "0.0 g".
The baby scale is ready for weighing.
4. Place the baby in the cradle.
The weight value will be flashing for 3 times and held. The screen displays the baby's stabilised weight.
Under the baby mode, if your baby is moving the value will be automatically hold.
5. Keep pressing  / TARE for 3 seconds to turn off the scale.
When it is not being used, the baby scale automatically switches off after approximately 20 to 30 seconds.

NORMAL MODE

Your Terrailon baby scale can also be used to weigh toddlers (up to 20 kilograms), by using only the glass platform without the baby tray.

1. Select first normal mode [OFFICE] with the mode switch on rear side of LCD panel.
2. Press  / TARE button to turn the power on, all digital will show on LCD for 5 secs. then "0" g shows on the LCD.
The scale is ready for weighting.
Simply have your child stand on the platform. The weight will be displayed.
3. Keep pressing  / TARE for 3 seconds to turn off the scale.
Or the power will be shut off automatically if there is not any button been touched in 4 to 8 minutes.

TERRAILLON EVOLUTIVE BABY SCALE

TARE FUNCTION

The tare function allows you to weigh your baby comfortably lying on a towel without including the weight of the towel.

- Switch on the scale by pressing the  / TARE button.
- Put the towel on the cradle. The weight value of the towel will be flashing for 3 times and held.
- Press the  / TARE button shortly, after 1 to 2 seconds the screen displays "0.0 g".
- Place your baby comfortably on the towel.

Only the baby's weight is displayed.

(When you remove both the baby and the towel, a negative weight may appear. This is the corresponding weight of the towel in negative.)

Notice:

1. Under  / TARE Weight function, the total weight cannot exceed the capacity of the scale.
2. If objects are taken out after  / TARE setting, the LCD will show the weight of the removed objects (in minus). Press  / TARE button to reset to zero.

ZERO CALIBRATION

Press  / TARE button to reset to zero.

WARNING INDICATORS

Problems	Reason	Solution
-----	The baby scale is overloaded (weight over 20 kg).	Remove your child immediately or you may damage the baby scale
LO	The battery is worn out.	Replace the battery immediately.

PROTECT THE ENVIRONNEMENT

Used batteries should be placed in the appropriate collection bin so they can be recycled. Do not mix different types of batteries, or use new and old batteries at the same time. At the end of its serviceable life, this product must be disposed of at a recycling centre for electrical and electronic equipment.

GUARANTEE

This product is guaranteed against defects in materials or manufacturing. During the guarantee period, any such defects will be repaired free of charge (proof of purchase must be shown in the case of a complaint under guarantee). The guarantee does not cover damage resulting from improper use or negligence. Claims during the first year of guarantee shall be addressed to the store where you purchased the product. From the second year of guarantee, claims should be addressed to Terrailon UK Ltd, 2 The Waterhouse, Waterhouse Street, Hemel Hempstead, Hertfordshire HP1 1ES. This commercial warranty is only applicable for products purchased from UK retail.

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Terraillon SAS France & Headquarters
Terraillon UK Ltd
Terraillon Asia Pacific Ltd
Terraillon Corp USA