
P0#1B2619 GBF-1266-C 说明书

印色: CMYK

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WEBCOACH style

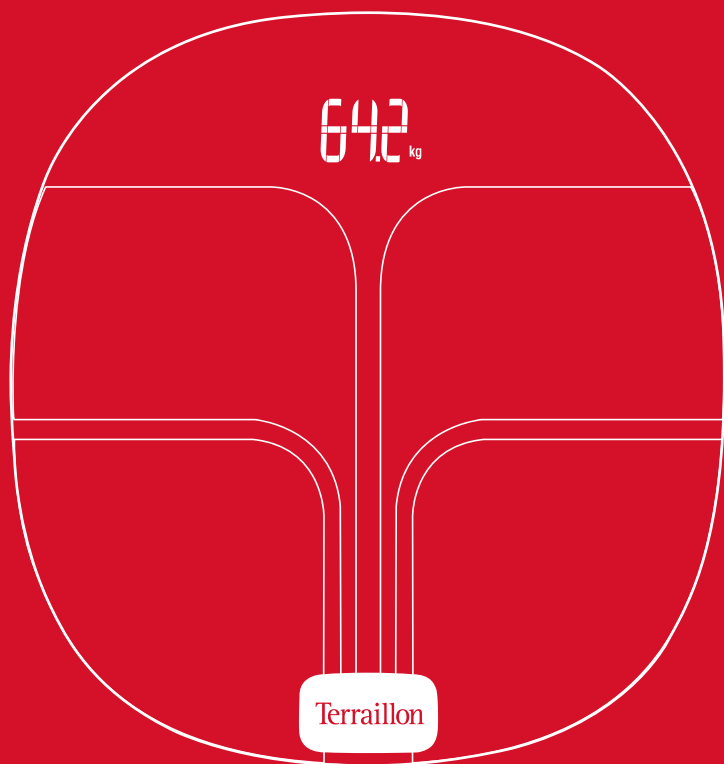USER MANUAL | GUIDE D'UTILISATION

NL Handleiding | **IT** Manuale d'uso | **ES** Manual del usuario

DE Bedienungsanleitung | **PT** Manual do utilizado

Terraillon®

www.terraillon.com

- EN** Thank you very much for selecting the Terraillon Glass Body Composition WEB COACH Style.
- FR** Nous vous remercions d'avoir choisi l'impédancemètre WEB COACH Style de Terraillon.
- NL** Bedankt dat u hebt gekozen voor de glazen lichaamsvetanalysator WEB COACH Style van Terraillon.
- IT** I nostri più vivi ringraziamenti per l'acquisto dell'analizzatore di massa grassa con piatto in vetro Terraillon WEB COACH Style.
- ES** Muchas gracias por elegir el analizador de masa grasa de cristal WEB COACH Style de Terraillon.
- DE** Vielen Dank, dass Sie sich für die Terraillon Körperfettanalysewaage WEB COACH Style entschieden haben.
- PT** Obrigado por escolher o Analisador de Massa Gorda de Vidro da Terraillon WEB COACH Style.

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Made for / Compatibilité

iPhone® 4S/5/5C/5S/6/6Plus
iPod® Touch 5th generation

iPad® 3/4/Air
iPad® Mini

Android 4.3
Bluetooth Smart 4.0 / Bluetooth Smart Ready

The range between the Analyzer and the Bluetooth end should be reasonably close, from 1 meter to 10 meters | **La distance** entre l'impédancemètre et le terminal Bluetooth doit être assez courte, entre 1 mètre et 10 mètres | **De afstand** tussen de analysator en de Bluetooth-ontvanger moet vrij klein zijn, van 1 tot 10 meter. | **La distanza** tra l'analizzatore e il dispositivo Bluetooth di destinazione deve essere ragionevolmente ridotta, da 1 a 10 metri | **La distancia** entre el analizador y el dispositivo Bluetooth debe ser razonablemente corta (de 1 a 10 metros) | **Der Abstand** zwischen Analysewaage und Bluetooth-Gerät sollte zwischen 1 Meter und 10 Meter betragen | **A distância** entre o Analisador e o terminal Bluetooth deve ser razoavelmente curta, entre 1 metro e 10 Metros

MAINTENANCE | MANUTENANCE

Onderhoud | Manutenzione | Mantenimiento | Wartung | Manutenção

Use water only. DON'T use propellant or other chemicals | **Utilisez** un chiffon humide uniquement. N'utilisez JAMAIS de détergents ou d'autres produits chimiques | **Gebruik** alleen water. GEEN drijfgas of andere chemische producten. | **Usare** con acqua. NON usare benzina o altre sostanze chimiche | **Limpie** con agua. NO use productos propelentes ni otro tipo de productos químicos | **Verwenden** Sie nur Wasser. Schmutz NICHT mit Treibmittel oder sonstigen chemischen Mitteln abwischen | **USE** água. NÃO gás ou outros produtos químicos

EN WARRANTY & ENVIRONMENT PROTECTION

This product is guaranteed for 2 years against material and manufacturing defects. During this period, such defects will be repaired free of charge (proof of purchase must be presented if the scale is under warranty). This warranty does not cover damage resulting from accidents, improper use or negligence. If you have a complaint please contact the store where you purchased the scale first.

Waste electrical products should not be disposed of with household waste. Please recycle where facilities exist. Check with your local authority or retailer for recycling advice.

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iPhone is registered trademark of Apple Inc.

Android is a trademark of Google Inc.

The Bluetooth® word mark and logos are registered trademarks owned by Bluetooth SIG, Inc.

Other trademarks and trade names are those of their respective owners.

Terraillon SAS France & Headquarters

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Terraillon Asia Pacific Ltd

4/F, Eastern Centre
1065 King's Road
Quarry Bay - Hong Kong
Tel: + 852 (0)2960 7200
customerservice@terraillon-asia.com

Terraillon Corp USA

contact@terraillon.fr

Terraillon®

www.terraillon.com

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Terraillon UK Ltd