QUICK START GUIDE

FR Guide de démarrage rapide
IT Guida rapida introduttiva

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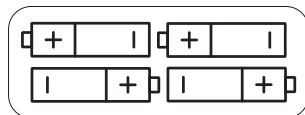
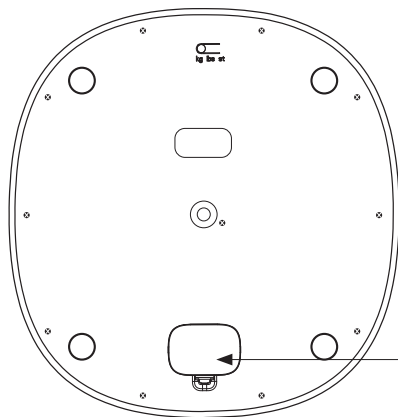
INSERT BATTERIES INTO THE DEDICATED COMPARTMENT.

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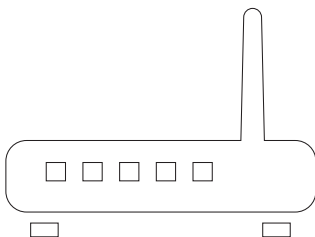
GO OVER TO YOUR INTERNET ROUTER AND MAKE SURE YOU KNOW YOUR WI-FI PASSWORD. THEN CONNECT YOUR SMARTPHONE TO YOUR WI-FI NETWORK.

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DOWNLOAD WELLNESS COACH – MY HEALTH MOBILE APPLICATION.

- FR** Téléchargez l'application mobile WELLNESS COACH – MY HEALTH.
NL Download de mobiele applicatie WELLNESS COACH – MY HEALTH.
ES Descarga la aplicación móvil WELLNESS COACH – MY HEALTH.
IT Scaricare l'app mobile WELLNESS COACH – MY HEALTH.
PT Descarregue a aplicação móvel WELLNESS COACH – MY HEALTH.
DE Laden Sie die mobile App WELLNESS COACH – MY HEALTH herunter.


CREATE YOUR ACCOUNT OR LOG IN WITH YOUR CREDENTIALS.

- FR** Créez votre compte ou connectez-vous avec vos identifiants.
NL Maak een account aan of meld u aan met uw inloggegevens.
ES Cree una cuenta o conéctese con sus credenciales.
IT Creare un account o accedere con i propri dati di accesso.
PT Crie a sua conta ou ligue-se com os seus identificadores.
DE Erstellen Sie Ihr Konto oder melden Sie sich mit Ihrem Login an.

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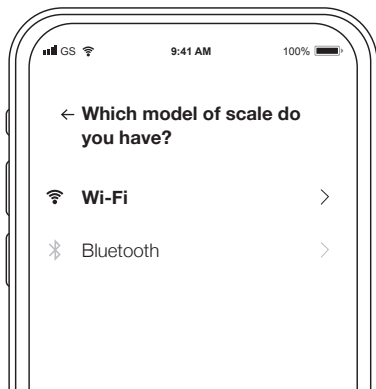DURING THE CREATION OF YOUR ACCOUNT, AND WHEN THE APPLICATION INVITES YOU TO DO SO, SELECT THE WI-FI BATHROOM SCALES TO START THE PAIRING PROCEDURE.

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NOTES: IF YOU ARE CONNECTED WITH A PRE-EXISTING ACCOUNT, PRESS THE BUTTON AS INDICATED, THEN "ASSOCIATE A NEW DEVICE". WHEN THE APPLICATION INVITES YOU TO DO SO, SELECT THE WI-FI BATHROOM SCALES TO START THE PAIRING PROCEDURE.

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
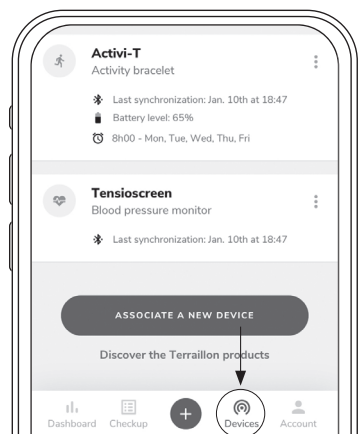
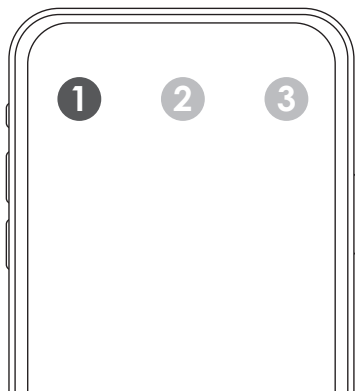
FOLLOW THE INSTRUCTIONS ON THE MOBILE APPLICATION.

- FR** Suivez les instructions dans l'application mobile.
- NL** Volg de instructies in de mobiele applicatie.
- ES** Siga las instrucciones de la aplicación móvil.
- IT** Seguire le istruzioni dell'applicazione mobile.
- PT** Siga as instruções da aplicação móvel.
- DE** Folgen sie den anweisungen in der mobilen anwendung.

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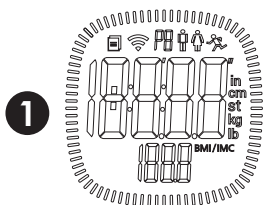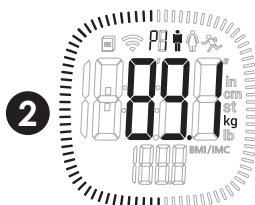
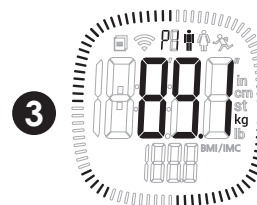
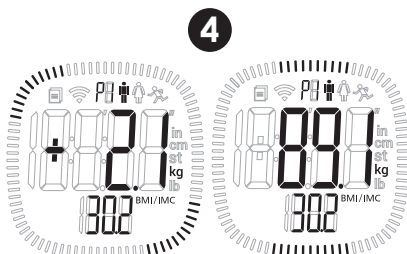
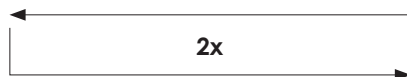
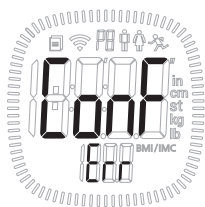
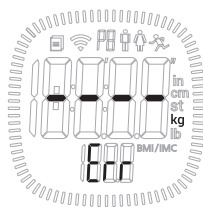
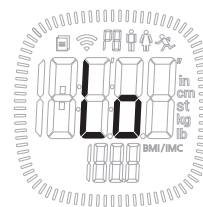
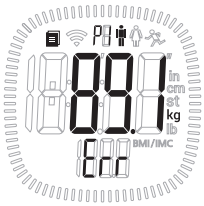
EN When the scale does not recognize user, two profiles are displayed (here P1 et P8).
 Press ▲ to select profile on the left, ▼ for the one on the right.

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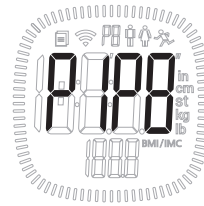
EN

Please read this section carefully to familiarize yourself with features and operations before using the unit. The warning signs and the sample icons shown here are listed in order for you to use this product safely and correctly as well as to prevent product damage, risk and injury to you or others. This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.

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EN Indoor use only, at temperatures of between 10°C and 40°C.

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EN Do not place the device in contact with liquids.

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EN

PRODUCT SPECIFICATION

Wi-Fi 802.11 b/g/n (2.4 Ghz compatible)
Weight & BMI (5 kg-180 kg)
8 user memories (with automatic recognition)
Battery powered (4 x AAA alkaline batteries (not included))
Dimensions: 32.4 x 34.2 x 2.3 cm

MAINTENANCE

Do not disassemble the scale; it contains no user serviceable parts. Do use a wet soft cloth, dipped into water and wrung out, to wipe the dirt. Then use a dry soft cloth to dry up the device. DON'T use propellant, abrasive or other chemicals.

PROTECT THE ENVIRONMENT

Dispose of used batteries in a designated container so that they can be collected and recycled. After life, dispose of this product at a designated waste recycling point.

WARNING

Do not expose batteries to high temperatures on account of the risk of liquefaction, vaporisation and explosion. Please ensure the batteries are replaced with batteries of the same type (otherwise your device may be damaged).

DECLARATION OF CONFORMITY

Hereby, Terrailon SAS declares that MASTER FORM complies with the essential health and safety conditions. This product meets the essential requirements and other relevant provisions of the directive RED 2014/53/UE. The complete declaration of conformity is available on <https://terrailonhelp.zendesk.com/hc/en-us> (Declaration of Conformity section).

WARRANTY

This product is guaranteed against defects in materials or manufacturing for 2 years. During this period, any such defects will be repaired free of charge (proof of purchase must be shown in the case of a complaint under warranty). The warranty does not cover damage resulting from improper use or negligence. In the event of a claim, first contact the shop where you purchased your product.

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FCC ID: 2AM6B-14711

FCC STATEMENT

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions:

1. This device may not cause harmful interference, and
2. This device must accept any interference received, including interference that may cause undesired operation.

CAUTION: The user is cautioned that changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

FCC RADIATION EXPOSURE STATEMENT

This equipment complies with FCC radiation exposure limits set forth for an uncontrolled environment.

This transmitter must not be co-located or operating in conjunction with any other antenna or transmitter

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

RF WARNING STATEMENT

The device has been evaluated to meet general RF exposure requirement, the device can be used in portable exposure condition without restriction.

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Terraillon SAS France & Headquarters

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16 Great Queen Street - Covent Garden
LONDON
WC2B 5AH
support.uk@terraillon.com

Terraillon Asia Pacific Ltd

4/F, Eastern Centre
1065 King's Road
Quarry Bay – Hong Kong
Tel: +852 (0)2960 7200
customerservice@terraillon-asia.com

Terraillon Corp USA

contact@terraillon.fr

MASTER FORM

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Terraillon UK Limited