PO#1SB002722A_GBF-1716-B1-004_PN15287_说明书

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
参考PO#1SB002722

BODY CONNECT

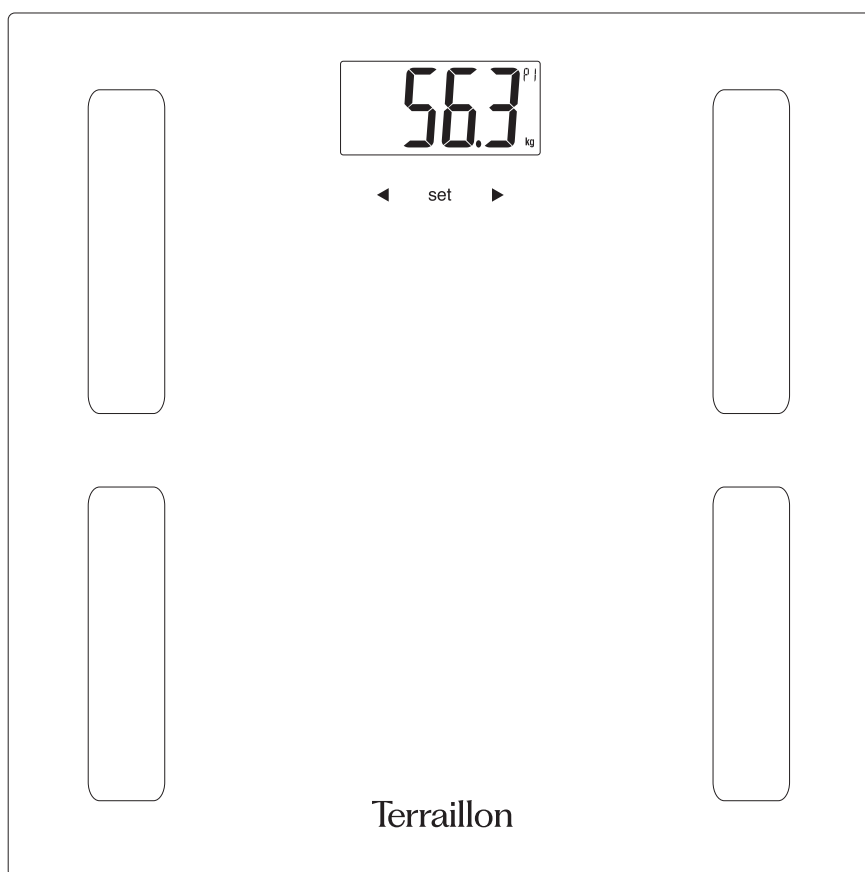USER MANUAL | GUIDE D'UTILISATION

Terraillon®

www.terraillon.com

Made for / Compatibilité

iOS 8.0 & +

Android 4.3 & +
Smartphone only / No tablet

iPhone is registered trademark of Apple Inc.

Android is a trademark of Google Inc.

The Bluetooth® word mark and logos are registered trademarks owned by Bluetooth SIG, Inc.

Other trademarks and trade names are those of their respective owners.

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
Thank you very much for selecting the Terraillon Glass Body Fat Analyzer BODY CONNECT.

Don't use this product if you wear a pacemaker or implantable electronic device • Any information provided by this device is in no way meant to treat, cure or prevent any disease or illness from happening.

PRODUCT OVERVIEW

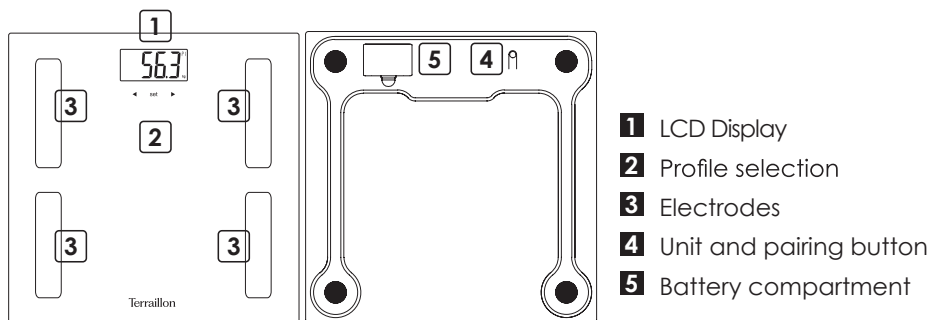FIRST USE

- Insert the batteries (3xAAA).
- You will need your smartphone to install your bathroom scales.
- Download the Wellness Coach - MyHealth app, then open it.

- Follow the instructions that appear on your smartphone. Ensure you enter all the information requested.
- Your smartphone must remain near to the bathroom scales throughout the installation process.
- It's important to bind your account with the scale for automatic recognition of user during weighing. You can pair the scale with 8 smartphones maximum. Before you first weighing, please select your profile number on the scale.

When the scale does not recognize user, two profiles are displayed (here P1 and P8). Press ◀ to select profile on the left, ▶ for the one on the right.

If you have any difficulties or need more information, you can find our help centre by scanning this QR code with your smartphone. If you cannot find the answer to your question, you can also contact our After Sales Service: support.uk@terraillon.com

START MEASURING

A

1. You must have bare feet to be able to analyse your body mass.
2. Open the MyHealth application.
3. Step onto your bathroom scale. Maintain contact with the electrodes until the following progress bar disappears from the screen.
 
4. If the application is not open, the data will be saved on the bathroom scale.
5. The symbol  will then be displayed on the screen.
6. To send the data, open the application on your smartphone, select your profile on the scale, then press SET. Otherwise, it will be sent the next time you weigh yourself if you follow step A.

Be sure to always weigh yourself on a hard and flat surface.

MAINTENANCE

Use water only. DON'T use propellant or other chemicals.

WARRANTY & ENVIRONMENT PROTECTION

This product is guaranteed for 2 years against material and manufacturing defects. During this period, such defects will be repaired free of charge (proof of purchase must be presented if the scale is under warranty). This warranty does not cover damage resulting from accidents, improper use or negligence. If you have a complaint please contact the store where you purchased the scale first.

Waste electrical products should not be disposed of with household waste. Please recycle where facilities exist. Check with your local authority or retailer for recycling advice.

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Terraillon SAS France & Headquarters
Terraillon Asia Pacific Ltd
Terraillon Corp USA