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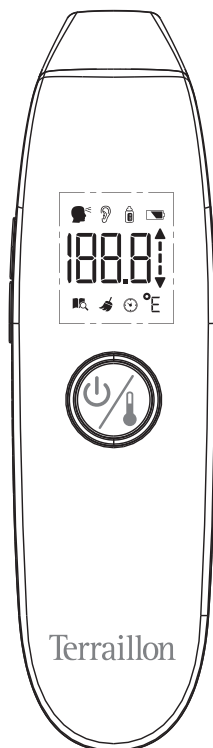
3 IN 1 AURICULAR THERMOMETER


Dear customer,

You have just bought this product and we thank you for your purchase. We hope you enjoy using the product. In order to be completely satisfied, we recommend that you read these directions carefully.

This thermometer is intended for home use by customer only. The measure site is ear canal.

MAINTENANCE

- Clean the thermometer with a soft, dry cloth.
- Clean the temperature probe with a cotton bud soaked in alcohol and dry immediately with a dry clean cotton bud. After cleaning, allow at least 20 minutes drying time before taking temperatures.
- Do not submerge the thermometer in water or any other liquid.
- Do not rinse the thermometer in water.
- Do not use water directly on the LCD screen.
- Store the thermometer in the dry, clean place.
- If the thermometer is not used for a long time, remove the batteries.

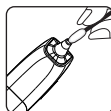SAFETY

- Keep out of reach of children.
- Do not use this appliance for any other purposes than those for which it was designed (human body temperature).
- Use of this thermometer is for monitoring purposes only and is not intended as a diagnostic tool. Make sure that any unusual readings are communicated to your physician. Under no circumstances should you change your prescribed medicine doses.
- Use of this thermometer is not intended as a substitution for consultation with your physician. Measurement results are for reference only. Contact your physician if you have or suspect any health concerns.

3 IN 1 AURICULAR THERMOMETER

RECOMMENDATION

- Store and use the thermometer in the same room.
- Do not use in areas where it may be subject to room temperature variations.
- Determine your normal temperature range by taking your temperatures when healthy. These readings can be used as a form of reference to compare with future readings when ill.
- Before taking a temperature, ensure that the forehead is free of hair and sweat, which may affect the reading. Also ensure that the probe is clean and undamaged.
- Earwax in the ear canal or lipid on the forehead may cause an inaccurate measurement, so make sure the ear canal and forehead is clean. The patient may be the proper operator of the thermometer. Only use the thermometer as described in this user manual. Any other use is deemed improper and may result in damage to property or even personal injury. When the device is used to measure the temperature of a child, it has to be operated by an adult. Adults can measure their own temperature.
- Do not take temperature within 30 minutes of eating, exercising, or taking a bath.
- Wait one minute between each temperature reading. Two consecutive readings may give slightly different measurements. Should this occur, take the average of the two.
- Do not attempt to disassemble the thermometer.
- Keep thermometer away from direct sunlight and store well away from any source of heat in a clean, dry location.
- Don't store the thermometer to extremes temperature or extremes humidity .Failing to do so may cause inaccuracy. Don't operate the thermometer under extremes temperature or extremes humidity. Failing to do so may cause inaccurate measurement.
- Do not drop or break the thermometer. Do not use if broken. Do not disassemble.
- Don't use a mobile phone nearby when the thermometer is working or do not use in an area with possible electromagnetic interferences.
- If these recommendations are not respected, the readings may be inaccurate.
- Temperature will vary depending on the part of the body. This is why you should not compare rectal and forehead temperature readings, for example.
- Temperature will also vary depending on the person, age, and time of day.
- Do not service and maintain the device while it is in use. Do not use any accessories which have not specified in this instruction manual. Keep the thermometer out of the reach of children and to avoid inhalation or swallowing of small parts. Do not allow children to take their temperatures unattended. If the temperature of the storage area differs greatly from that of the measuring area, please wait the thermometer temperature to equalize to the room temperature about 30 minutes before measuring. Please avoid mechanical shock, otherwise the thermometer will be damaged.

3 IN 1 AURICULAR THERMOMETER

BODY TEMPERATURE

Normal body temperature is a range. A person's normal temperature range tends to decrease with age. The following table shows normal temperature ranges by age:

0 – 2 years	36.4 – 38 °C	97.5 – 100.4 °F
3 – 10 years	36.1 – 37.8 °C	97.0 – 100.0 °F
11 – 65 years	35.9 – 37.6 °C	96.6 – 99.7 °F
> 65 years	35.8 – 37.5 °C	96.4 – 99.5 °F


The range of normal temperature varies from person to person and can be influenced by many factors such as time of day, measure from different sites, level of activity, medications, emotion and so on. So we recommend that you practice with the Thermometer on yourself and family members when you are healthy. This way you know how the thermometer works and can feel more confident of the measurements you take when a family member is ill.

PACKAGE COMPONENTS


NO.	DESCRIPTION	QUANTITY
1	Ear and Forehead Thermometer	1 piece
2	1.5V battery AAA	2 pieces
3	User manual	1 piece

3 IN 1 AURICULAR THERMOMETER

PRODUCT DESCRIPTION

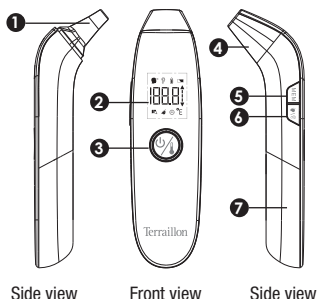
- ① Infrared sensor
- ② LCD screen
- ③ Start/stop button
- Taking a temperature
- ④ Probe cap
- ⑤ 20-reading recall button
- ⑥ Temperature mode button
- ⑦ Battery compartment

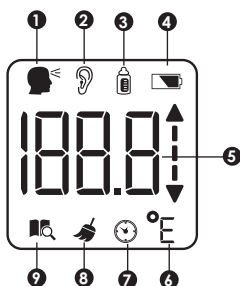- ① Forehead temperature mode
- ② Ear temperature mode
- ③ Object temperature mode
- ④ Battery empty symbol
- ⑤ Temperature on LCD screen
- ⑥ Degrees °C, °F, E
- ⑦ Temperature taking recall
- ⑧ Clean probe symbol
- ⑨ Memory symbol

INSERTING & CHANGING BATTERIES

- To insert or change the 2 LR3 batteries (supplied with the thermometer), open the battery compartment by sliding off its cover, insert the batteries ensuring that the negative and positive ends are in the right place, and replace the cover before use.
- When the low battery indicator  appears on the LCD, or the thermometer does not function at all, you should replace new battery as soon as possible.

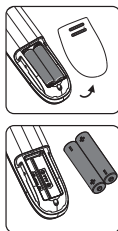3 IN 1 AURICULAR THERMOMETER

CHOOSING THE TEMPERATURE MEASURING UNIT (C° OR F°)

- Press the on/off button to turn on the thermometer. You can switch between °C (Celsius) and °F (Fahrenheit) by holding down the [MEM] button.

USING THE PRODUCT

FOREHEAD TEMPERATURE

Check that the cap is in place for the front temperature measurement.

- Press the mode button to select 'take forehead temperature'.
(Check that the probe is clean before taking a temperature).
- Place the probe on the forehead.
- Press the 'take temperature' button for 2 seconds.
- Keep the thermometer on your forehead until two beeps.
- Then remove the thermometer from the forehead.
- The temperature will be displayed and the colour red or green depending on the result.
- After 6 seconds, you can retake your temperature.
- Temperature is automatically memorised.


The forehead temperature may be affected by external climatic conditions so, in case of doubt, do not hesitate to take the ear temperature.

Red: 38.0°C (100.4°F) ~ 42.9°C (109.3°F)

Green: 34.0°C (93.2°F) ~ 37.9°C (100.3°F)

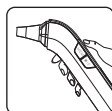CAUTION: External factors may influence forehead temperature particularly when an individual has:

- been exposed to very hot or very cold temperatures.
- been recently swimming or bathing.
- had their forehead covered.

In these cases, let the individual stay away from the above scenario and station in a room for 30 minutes prior to taking a temperature.


3 IN 1 AURICULAR THERMOMETER

TAKING A TEMPERATURE IN THE EAR

- Remove the probe cap.
- Press the mode button to select 'take ear temperature'. (Check that the probe is clean before taking a temperature).
- Place the probe in the ear canal.
- Press the 'take temperature' button for 2 seconds.
- Keep the thermometer in your ear canal until you hear two beeps.
- Then remove the thermometer from the ear canal.
- The temperature will be displayed and the colour red or green depending on the result.


After 6 seconds, you can retake your temperature.

If you take the temperature in your right ear, the reading may be different from a temperature taken in your left ear. Therefore, you should always take your temperature in the same ear.

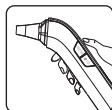CAUTION:

In order to obtain a precise measurement, ensure the correct positioning of the probe in the ear so that the probe tip senses the heat from the eardrum. In order to do this, you can pull the ear to straighten the ear canal and slightly rotate the probe inside. (It is not recommended to use the ear function of the thermometer on a baby under 12 months old.)

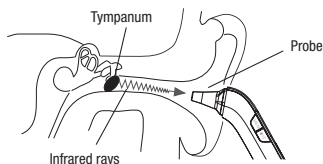Easy interpretation of the results

In both modes, if the temperature is lower or equal to 37.9° C, the green light comes on for around 5 seconds. If the temperature is higher than 37.9° C, the red light comes on.

- Temperature higher than 37.9° C - high temperature: the red light comes on.
- Temperature lower or equal to 37.9° C - normal temperature: the green light comes on.

3 IN 1 AURICULAR THERMOMETER

TAKING THE TEMPERATURE OF AN OBJECT




- Remove the probe cap.
- Press the mode button to select 'take temperature of an object'.
- Always place the probe 1 to 5 cm from the object.
- Press the 'take temperature' button and hold it down for 2 seconds.
- Then remove the thermometer from the object.
- The temperature reading will be displayed.

CAUTION:

Make sure don't use this mode to measure body temperature.

CAUTION:

With object measurement, the current surface temperature of the object is displayed. It can be different from its internal temperature, especially if the surface is exposed to direct sunlight or a draught.

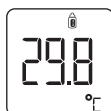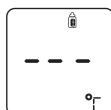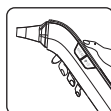MEMORY

- There are a total of 20 memory data for recording forehead temperature measurements.
- The current measurement is always stored in the last storage space. When all storage spaces have been occupied, the oldest measurement is always deleted from the memory.





To display the last memorised result, press the memory recall button «MEM». The last temperature reading will appear along with the memory number.

To display the previous memorised results, press the memory recall button «MEM» again and the memory numbers will be displayed.

Then press the "Mode" button to return to temperature-taking mode.

Compare against your last reading.

The curve lets you see quickly if your temperature is higher or lower than your last reading.

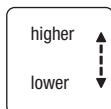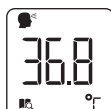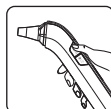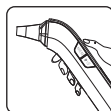Terraillon®

3 IN 1 AURICULAR THERMOMETER

CLEANING TIP

For reminding user to clean the probe tip, this thermometer has a special design --- cleaning tip.

When you have not used the thermometer for a long time or done measurement after multiple times, the clean symbol appears «» to remind the user to clean the infrared sensor. The thermometer will return back to original mode after any button is pressed, and the backlight stop flashing.

TROUBLESHOOTINGS

Troubles	Checklists	Counter measures
	The measured temperature is lower than 32.0°C/89.6°F or higher than 42.9°C/109.3°F Please check the operation method	Follow User's Manual for proper measurement.
	Hardware problem	Contact the retailer or customer services
	Operating temperature is out of the range	Use the thermometer in the range of operating conditions.
The measurement is not accurate or if there is any doubt on the measured result	Please check if the infrared sensor is clean or not	Clean the infrared sensor with cotton swab according to user manual
	Please check the measuring way is correct or not	Ensure you have read the manual and know how to use the thermometer properly. Please keep the thermometer and patient in the measuring room at least 30 minutes before using Please take the measurement indoor Do not held the thermometer in your hand too long. It could affect the accuracy
	Please check if there is low battery symbol on LCD	Change new battery

3 IN 1 AURICULAR THERMOMETER

PROTECT THE ENVIRONMENT

Dispose of used batteries in a designed container so that they can be collected and recycled. Do not mix different types of batteries. Do not mix new and used batteries. After life, dispose of this product at a designated waste recycling point.

GUARANTEE

This product is guaranteed against defects in materials or manufacturing. During the guarantee period, any such defects will be replaced free of charge (proof of purchase must be shown in the case of a complaint under guarantee). The guarantee does not cover damage resulting from improper use or negligence. In the event of a claim, first contact the shop where you purchased the product.

Applied Standard.

- IEC/EN 60601-1 (Safety)
- IEC/EN 60601-1-2 (EMC)

Keep these instructions in a safe place for future consultation.

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SPECIFICATIONS

Power supply:	DC 3V (2 x AAA battery)
Body measuring range:	32.0°C ~ 43.0°C (89.6°F ~ 109.4°F)
Laboratory accuracy:	±0.2°C(±0.4°F) for 35.0°C (95.0°F) ~ 42.0°C (107.6°F) ±0.3°C (±0.5°F) for other range
Object measuring range:	20°C~60°C (68°F~140.0°F)
Object measuring accuracy:	±2.0°C(±4.0°F)
Displayed room temperature range:	5.0°C ~ 59.9°C (41.0°F~139.8°F)
Room temperature accuracy:	±1°C (±2°F)
Measure time:	1 second
Measure interval:	1 second
Dimension:	145.4(L)x40.0(W)x42.9(H)mm(±2mm)
Weight:	About 65.0g (without battery)
Forehead operating condition:	15.0°C~40.0°C (59.0°F~104.0°F) with a relative humidity of 15%~85% Atmosphere pressure: 700hPa-1060hPa
Ear and object operating condition :	10.0°C~40.0°C (50.0°F~104.0°F) with a relative humidity of 15%~85%, Atmosphere pressure: 700hPa-1060hPa
Storage and transport condition:	-20°C~55°C (-4°F~131°F) with a relative humidity of 15%~90% Atmosphere pressure: 700hPa-1060hPa

Reference body site: core

Operating mode: adjusted mode

Service life: 2 years

Service life of the battery: With a new battery (Carbon) approx. 10 months based on the use frequency at 5 times/day

Adjusted mode	Item	Group A1	Group A2	Group B	Group C
Clinical bias-1(ear)	Δ_{cp-1}	-0.275	-0.135	-0.121	-0.505
Clinical bias-2(forehead)	Δ_{cp-2}	-0.025	-0.012	0.067	-0.298
Uncertainty-1(ear)	L_A-1	0.999	1.541	1.436	1.375
Uncertainty-2(forehead)	L_A-2	1.358	1.461	1.302	1.160
Clinical repeatability(ear)	σ_r-1	0.160	0.176	0.190	0.225
Clinical repeatability (forehead)	σ_r-2	0.222	0.210	0.214	0.239

Caution:

This thermometer has been thoroughly tested and inspected to assure proper performance and operation!

Terrailon®

NORMALIZED SYMBOLS

	Applied part of type BF (The applied part is the probe tip)
	Disposal in accordance with Directive 2012/19/EU (WEEE)
	Complies with the European Medical Device Directive (93/42/EEC) and amended by directive 2007/47/EC requirements. Notified body TÜV Rheinland (CE0197)
	The name and the address of the manufacturer
	The name and the address of the Authorized EC-representative in Europe
	Caution
	Refer to Instruction Manual
IP22	The first number 2: Protected against solid foreign objects of 12,5 mm Φ and greater. The second number: Protected against vertically falling water drops when enclosure titled up to 15°. Vertically falling drops shall have no harmful effects when the enclosure is titled at any angle up to 15°, on either side of the vertical.
	Date of Manufacture

Electromagnetic Compatibility (EMC)

- This device should not be used adjacent to or stacked with other equipment and that if adjacent or stacked use is necessary, this device should be observed to verify normal operation in the configuration in which it will be used.
- Use of accessories other than those specified or provided by the manufacturer of this device could result in increased electromagnetic emissions or decreased electromagnetic immunity of this equipment and result in improper operation.
- Portable RF communications equipment (including peripherals such as antenna cables and external antennas) should be used no closer than 30 cm (12 inches) to any part of the device, including cables specified by the manufacturer. Otherwise, degradation of the performance of this equipment could result.
- When the operating environment is relatively dry, strong electromagnetic interference usually occurs. At this time, the device may be affected as follows:
 - the device stops output;
 - the device turns off;
 - the device restarts;

The above phenomenon does not affect the basic safety and essential performance of the device, and the user can use it according to the instruction. If you want to avoid the above phenomenon, please use it according to the environment specified in the manual.

Table 1

Declaration - electromagnetic emission		
The device is intended for use in the electromagnetic environment specified below. The customer or the user of device should assure that it is used in such an environment.		
Emissions test	Compliance	Electromagnetic environment - guidance
RF emissions CISPR 11	Group 1	The device uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
RF emissions CISPR 11	Class B	The device is suitable for domestic establishment and in establishment directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.
Harmonic emissions IEC 61000-3-2	Not applicable	
Voltage fluctuations/ flicker emissions IEC 61000-3-3	Not applicable	

Table 2

Declaration - electromagnetic emission			
The device is intended for use in the electromagnetic environment specified below. The customer or the user of device should assure that it is used in such an environment.			
Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment - guidance
Electrostatic discharge (ESD) IEC 61000-4-2	±8 kV contact ±2 kV, ±4 kV, ±8 kV, ±15 kV air	±8 kV contact ±2 kV, ±4 kV, ±8 kV, ±15 kV air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30 %.
Electrical fast transient/burst IEC 61000-4-4	± 2 kV for power supply lines ± 1 kV for input/output lines	Not applicable	Mains power quality should be that of a typical commercial or hospital environment.
Surge IEC 61000-4-5	± 0.5kV, ± 1 kV line(s) to lines ± 0.5kV, ± 1 kV, ± 2 kV line(s) to earth	Not applicable	Mains power quality should be that of a typical commercial or hospital environment.
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	0 % UT; 0.5 cycle At 0°, 45°, 90°, 135°, 180°, 225°, 270° and 315° 0 % UT; 1 cycle and 70 % UT; 25/30 cycles Single phase: at 0° 0 % UT; 250/300 cycles	Not applicable	Mains power quality should be that of a typical commercial or hospital environment. If the user of the device requires continued operation during power mains interruptions, it is recommended that the device be powered from an uninterruptible power supply or a battery.
Power frequency (50/60 Hz) magnetic field IEC 61000-4-8	30 A/m	30 A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.
NOTE: UT is the a.c. mains voltage prior to application of the test level.			

Table 3

Declaration - electromagnetic immunity			
The device is intended for use in the electromagnetic environment specified below. The customer or the user of device should assure that it is used in such an environment.			
Immunity test	IEC 60601test level	Compliance level	Electromagnetic environment - guidance
Conducted RF IEC 61000-4-6	3V 0.15 MHz to 80MHz 6 V in ISM and amateur radio bands between 0.15 MHz and 80 MHz	Not applicable	Portable and mobile RF communications equipment should be used no closer to any part of device, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter. Recommended separation distance $d = 1.2 \sqrt{P} \text{ 150 KHz to 80 MHz}$ $d = 1.2 \sqrt{P} \text{ 80 MHz to 800 MHz}$ $d = 2.3 \sqrt{P} \text{ 80 MHz to 2.7 GHz}$ where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in meters (m). Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey,^a should be less than the compliance level in each frequency range. Interference may occur in the vicinity of equipment marked with the following symbol: 
Radiated RF IEC 61000-4-3	10V/m 80 MHz to 2.7 GHz	10V/m	

NOTE 1 At 80 MHz and 800 MHz, the higher frequency range applies.

NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

^a Field strengths from fixed RF transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which device is used exceeds the applicable RF compliance level above, device should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as re-orienting or relocating device.

^b Over the frequency range 0.15 MHz to 80 MHz, field strengths should be less than 3 V/m.

Table 4

Recommended separation distances between portable and mobile RF communications equipment and device			
The device is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of device can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and device, as recommended below, according to the maximum output power of the communications equipment.			
Rated maximum output power of transmitter W	Separation distance according to frequency of transmitter m		
	0.15 MHz to 80 MHz $d = 1.2 \sqrt{P}$	80 MHz to 800 MHz $d = 1.2 \sqrt{P}$	80 MHz to 2.7 GHz $d = 2.3 \sqrt{P}$
0.01	0.12	0.12	0.23
0.1	0.38	0.38	0.73
1	1.2	1.2	2.3
10	3.8	3.8	7.3
100	12	12	23
For transmitters rated at a maximum output power not listed above, the recommended separation distance d in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer NOTE 1 At 80 MHz and 800 MHz, the higher frequency range applies. NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.			

Model: TH1110A / THERMO FAST
Shenzhen Dongdixin Technology Co., LTD.
Floor 1-2, No.3 Building, Fanshen Xusheng Industrial Estate Xilixiaobaimang
518108 Nanshan District, Shenzhen P. R. China

Shanghai International Holding Corp. GmbH (Europe)
Eiffestraße 80, 20537 Hamburg Germany

Version: V1.0

Terraillon SAS France & Headquarters

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16 Great Queen Street - Covent Garden
LONDON
WC2B 5AH
support.uk@terraillon.com

Terraillon Asia Pacific Ltd

4/F, Eastern Centre
1065 King's Road
Quarry Bay – Hong Kong
Tel: +852 (0)2960 7200
customerservice@terraillon-asia.com

Terraillon Corp USA

contact@terraillon.fr

Terraillon®

www.terraillon.com

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Terraillon UK Limited